




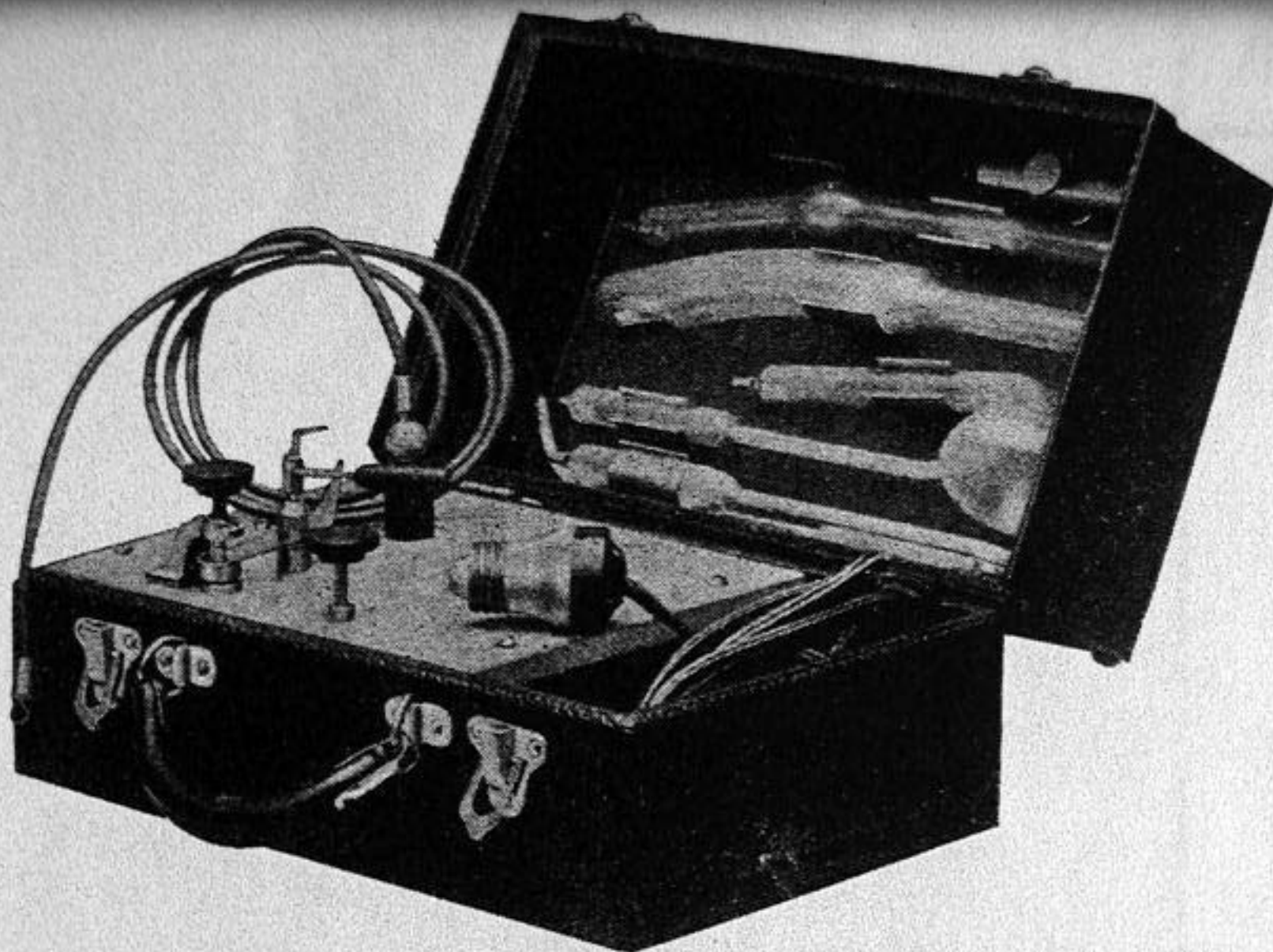

Type "F"

Complete Physician's Portable
High Frequency and
Violet Ray Generator.

Manufactured by

H. G. FISCHER

2009-11 Raymond St.

CHICAGO ILL.

There is no quackery or uncertainty about it, and it is not put forward as a "cure-all", but being a thorough cell massage and peculiarly healthy stimulant, its intelligent use is a wonderful aid in restoring the normal equilibrium wherever it has been disturbed and the physical condition impaired from almost any cause, and used as directed with ordinary discretion, its results are almost universally beneficial and never harmful.

The practitioner once initiated into the results obtained from the use of high frequency high potential violet ray currents such as generated by Type "F" finds their use absolutely indispensable.

We pride ourselves in the extreme high frequency produced by Type "F" Violet Ray Generator. We have a current here that the Physician can use with absolute safety to himself and his patient, and rest assured that he has done all that can be done in the shortest possible time.

Type "F" is not a toy and we want to call your attention to the fact that when Button E is pressed to bottom, you have a powerful current to deal with. A current that may aggravate or even burn your patient, when used where not indicated. (We refer chiefly to internal treatments.) Therefore we say play safe, doctor, and remove Button E entirely when you are doing delicate work. This will positively assure you that you are using currents of highest possible frequency.

No other outfit has this feature and the physician using other outfits is compelled to regulate as best he can with a current of medium power, which is neither high nor low enough in frequency to cover all cases.

Bear in mind that all internal treatments must be made with extreme high frequency to avoid discomfort to the patient and Plug "E" must be pulled out.

Another feature combined in every Type "F" is a Tesla Coil protector, which insures long life and enables us to guarantee it for twice the length of time of any other coil. Hence our statement it will outlive you.

Type "F" outfit is substantially built upon **White Italian Marble**, and is inclosed in a **plush lined** case with **steel Valence**. The outer trimmings being of **brass** heavily nickel plated.

GENERAL INSTRUCTIONS

For the use of high frequency violet rays generated by Type "F" outfit

We recommend the more detailed method of instruction as given in the following text-books.

High Frequency Manual.....	Eberhart...	\$2.50
High Frequency Currents.....	Strong.....	3.00
Currents of High Potential....	Snow.....	3.50

The following are a few suggestions for treating diseases taken from cases where physicians have obtained excellent results. Not all of them through high frequency alone, but as a valuable adjunct in the treatment thereof.

Abscesses Electrode No. 1. Keep in contact with skin and treat for 5 to 7 minutes with medium spark high frequency current.

Acne The fulguration electrode No. 20 is used for quickest results. Cauterize the larger pustules and treat the entire surface affected with the No. 1 Electrode for a few minutes.

Acne Rosacea Electrode No. 1 or No. 2. Medium frequency. Keep electrode away from surface to get cauterizing effect of spark. Duration of treatment to be decided by severness of case. Where much surface is to be covered lengthen time to 10 minutes.

Read Eberhart's Manual, pages 164 to 168.

Alopecia Electrode No. 1 or No. 11. Use medium frequency over the scalp for 5 minutes. If current is too strong use higher frequency and lengthen the time.

Anus Fissure of Electrode No. 2 is to be used. Treat about 4 to 6 minutes with high frequency current as strong as patient can bear without discomfort.

Diagnosis like Criticism is of little avail unless it travels hand in hand with Remedy, Treatment and Cure.

Headache When not caused by stomach apply Electrode No. 1 with medium spark high frequency over forehead and back over the ears to back of neck. Occasionally down the spine.

Hemorrhoids Use rectal electrode No. 2 or 12. Use delicate High frequency current. To remove use fulguration Electrode No. 20.

Hypertension Auto condensation produces most rapid results.

See Arteriosclerosis.

Hypotension Apply Electrode No. 1 to spine with heavy spark.

Impotence Use insulated rectal Electrode No. 12 across the genitals with medium to high frequency current, for 4 to 6 minutes.

Influenza Breathing through ozone generator brings immediate relief. Follow with treatment over sides of nose and in some cases use nasal electrode internally.

Lumbago Electrode No. 1 with sharp spark over and beyond painful area. Treat every other day for about 6 minutes.

Massage Electrode No. 1. Medium or low frequency as preferred by patient. The lower the frequency and the sharper the spark the shorter the time required for treatment. Through heavy gauze is ideal.

Moth Patches or Moles Fulguration Electrode No. 20 with medium or high frequency, Do not overtreat as 1 to 2 minutes suffices over a small area.

Nervousness Treat over spine with medium to high frequency. Daily treatments at first, followed by treatments every three days for about one month.

Neuritis Surface Electrode No. 1, heavy current high frequency. Relief almost instantly. Do not treat over same surface over 5 minutes, but treat twice daily if at all possible. Results are most pleasing even in chronic cases.

Paralysis Electrode No. 1 with sharp spark of medium frequency. Apply over muscles in contact and by patting. Raising the Electrode about 1-4 to 1-2 inch above the surface. Treatments should be made daily and be of 4 to 8 minutes duration.

Prostatic Diseases Special Electrodes are made for treating diseases of the prostate, but excellent results have been obtained by use of Electrode No. 12, in reducing enlarged prostate. A mild medium frequency is used for 5 to 10 minutes.

Pyorrhea Alveolaris Apply Special Dental Electrode No. 19 to gums for 3 to 5 minutes, using high frequency and in direct contact with the gums.

Rheumatism Combine treatments as given in Neuritis and Paralysis.

Rodent Ulcer Best results are obtained with X-Ray. Fulguration Point No. 20 has been used to destroy the affected tissue.

Skin Diseases The Dermatologist will find many uses for the high frequency current due to its cauterizing, escharotic and antiseptic properties. Read any good book on the subject of high frequency or high potential currents.

Sciatica Body Electrode Medium frequency. Apply for 5 to 7 minutes.

Stiff Joints Apply as for Neuritis with longer treatment if deemed necessary.

Throat Diseases Use Throat Electrode No. 5 or No. 18 with mild current high frequency. Finish by treating above and below larynx with Electrode No. 33. Duration of treatment about 8 minutes.

Torticollis (Stiff Neck) Treat back of ear and down back and sides of neck. Use Electrode No. 1 with sharp spark medium frequency; patting the surface gently. Give treatments of 10 minutes or longer.

Trachoma Apply eye electrode No. 8 with highest frequency. Treat over closed eye lid for one to 3 minutes.

Ulcers (Chronic) Treat through gauze with as sharp a spark as can be tolerated, using Electrode No. 2 or surface No. 1. Ozone generator can be used to advantage by employing the spray.

Uterine Diseases Use Electrode No. 13 and apply with highest frequency for short periods. Beneficient results are always obtained.

Don't let your patient sit close to a radiator while treating.
 Don't let your electric light burn too near your patient. Remember you charge him with 60,000 volts and the discharge of current from his body to the surrounding air will soon burn out your lamp.
 Don't expect results too soon.
 Don't limit the use of your coil to one or two cases only, write for latest text books.
 Don't blame your coil if you cannot get a spark, look at your own wiring.
 The coil is made under a guarantee and is carefully inspected.
 And don't forget we have other electrical goods of equal quality.

PRICE LIST

Type "F" outfit complete.....	\$ 30.00
Type "F" outfit with addition of Cautey circuit.....	32.50

VACUUM ELECTRODES

No. 1. Surface.....	\$.75
2. Rectal.....	.75
3. Vaginal.....	.75
4. Nasal.....	.75
5. Throat.....	.75
6. Urethral.....	.75
7. Eye, Single.....	.75
8. Eye, Double.....	1.00
9. Metal disc body electrode.....	3.00
10. Dental, surface, insulated.....	1.50
11. Comb electrode.....	1.50
12. Rectal, insulated.....(Prostatic).....	1.50
13. Vaginal, insulated.....	1.50
14. Nasal, insulated.....	1.50
15. Post nasal, insulated.....	1.50
16. Urethral, insulated.....	1.50
17. Ear, insulated.....	1.50
18. Tongue, insulated.....	1.50
19. Dental, insulated.....	1.50
20. Fulguration.....	2.00
21. Ozone generator.....	7.50
22. Ozone generator complete with bulb.....	10.00
30. Corn Electrode.....	.75
33. Throat Electrode (outer).....	.75
50. Universal handle insulated.....	1.00
51. High frequency cord (heavy).....	.75
52. Metal handle for indirect treatments, extra long.....	.65
53. Glass handle for indirect treatments.....	1.00
301. Cautey handle adjustable to 5 positions.....	1.50
306. Cautey Cords per pair.....	1.00
305. Cautey knives any shape.....	.75
Complete set of SNOW'S electrodes Nos. 1, 2, 3, 4, 5, and 6.....	3.50

We also Manufacture and carry in stock a complete line of Electro Therapeutic Apparatus, which includes the following;

Auto Condensation Treatment Outfits

Auto Condensation Pads

Diathermic Treatment Outfits

X-Ray Plates and Equipment

High Frequency Outfits

X-Ray Generators

Thompson-Plaster Electrical Cabinets

Diagnostic Light Regulators

Portable Combination Outfits

Cautery Transformers

Cautery Equipment

Therapeutic Lamps




Head Lamps

Dental High Frequency Outfits

Dental Wall Plates

Ask your nearest dealer. If he will not supply you, send your order direct to us.

50

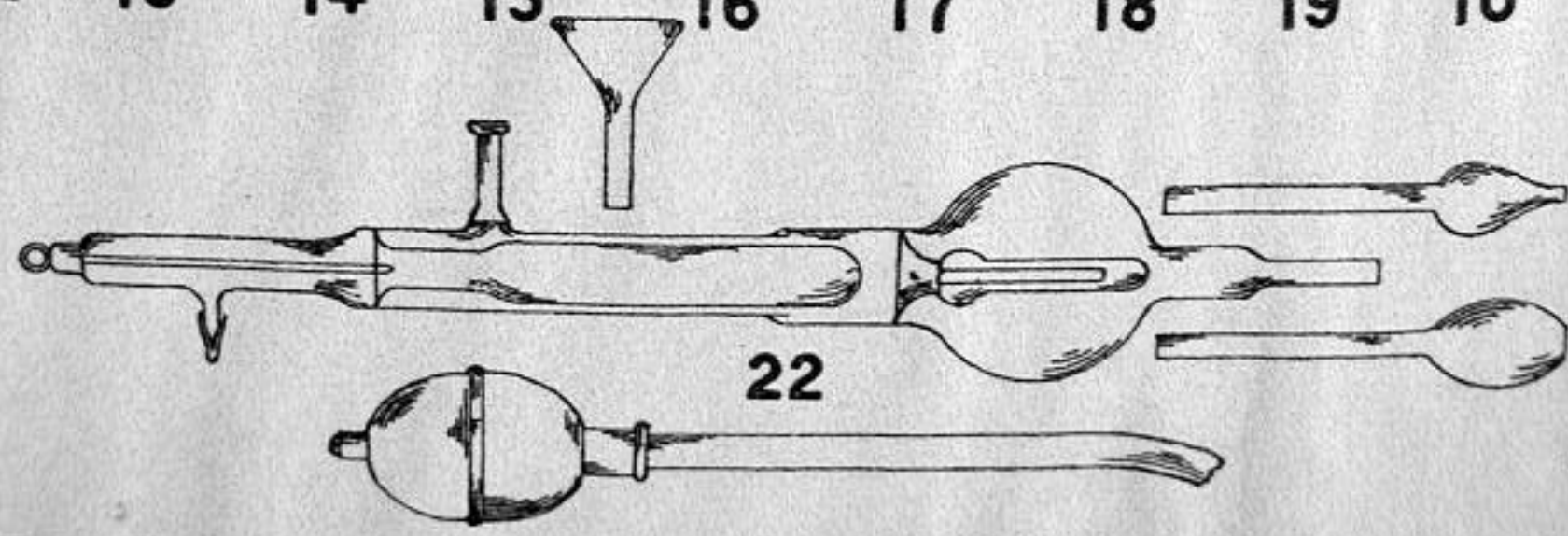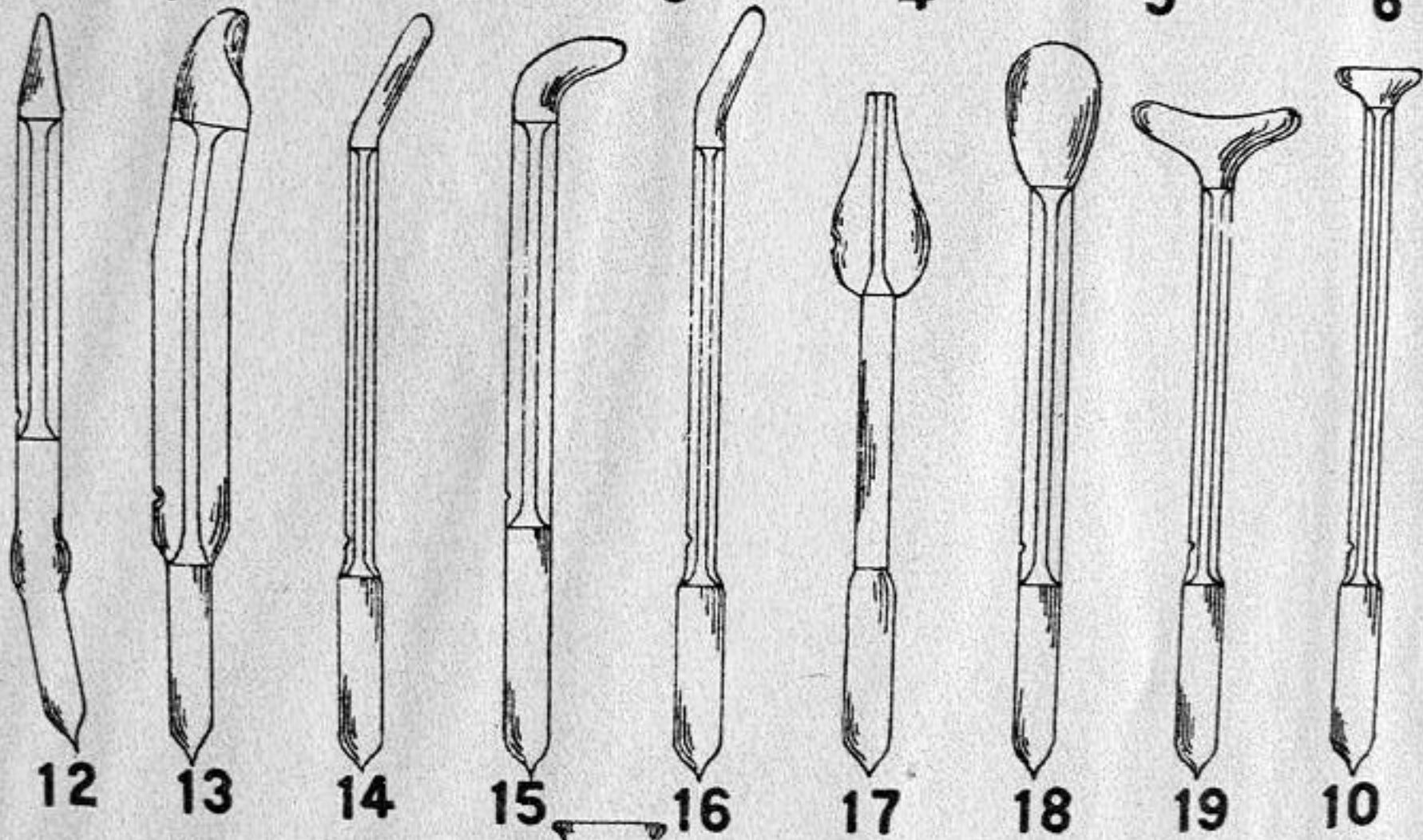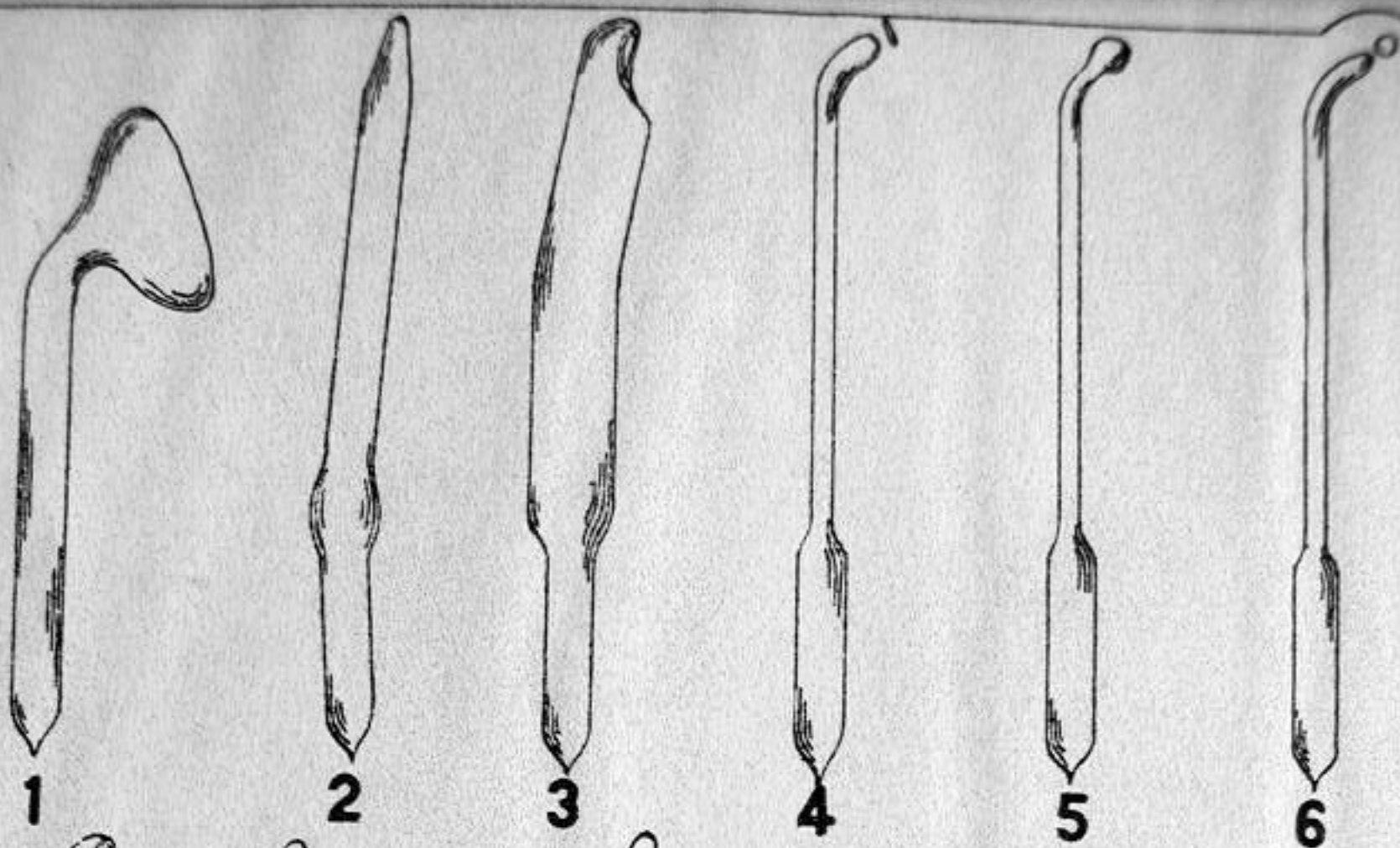

H



TYPE "F" COIL

is guaranteed for a
period of 2 years
from date of
purchase
against
inherent
defects.

Whether Intelligence rules
all things or Chance rules,
Man must use the Power
that he has and
be Tranquil.